


# Spring2

## SPRAY-LIGHT MICROMOTOR FOR PODIATRY AND DENTAL LAB



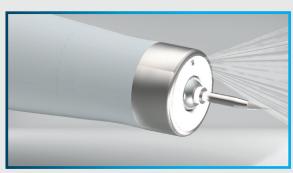

- Integrated SPRAY system with optimally nebulized jet and accurate adjustment.
- > High luminosity LED light.
- Control unit with water tank, easy and practical control-panel with knobs for spray and speed adjustment.
- > Handy and light handpiece with fast chuck changing.
- > Brushless micromotor of the newest generation, high torque and speed performance.

Detachable handpiece from the tube for service purpose



#### **Specifications:**

Adjustable speed up to 40.000 rpm

Adjustable spray jet, flow rate from 0 to 7 ml/min

Reverse rotation function

Automatic spray and/or light function at the motor launching

Automatic tube filling

Specifications LED: 130 lumen, 5000 °K

Control-unit: dimension (WxDxH) 235x311x118 mm (tank included), weight 3,7 kg

Tank capacity: 280 ml

Handpiece: length 14 cm, Ø min 2 cm, Ø max 2,5 cm, weight 170 g

On/off motor pedal







The High Quality Made in Italy



### **MARIOTTI & C. SRL**

#### Micromotors manufacturer since 1973

Via Seganti 73 - 47121 Forlì - Italy Tel +39 0543 474105 - Fax +39 0543 781811 info@mariotti-italy.com www.mariotti-italy.com







